

Pi - 102

The affordable multifunctional luminometer

Luminometer



Key features

- Portable luminometer with interchangeable counting chamber
- Suitable for chemi- and bioluminescence
- Sensitivity down to 0.5 fmole ATP
- Software for data handling included
- Programmable cut-off limits
- Battery/Mains operated
- RS232 serial port for PC/Printer
- Robust

User benefits

- Suitable for cuvettes and all-in-one reagent devices
- Multiple applications
 - a) ATP bioluminescence
 - surface contamination e.g. hygiene monitoring
 - aqueous contamination e.g. process waste water
 - b) Chemiluminescence
 - e.g. commercial and R+D diagnostic tests
- Most sensitive portable luminometer available
- Flexible data management options
- Objective automated interpretation of results
- Independent of location
- Stand-alone or PC compatible
- High quality and reliable



hygiena
INTERNATIONAL


With cuvettes



Quick and easy change-over



With devices

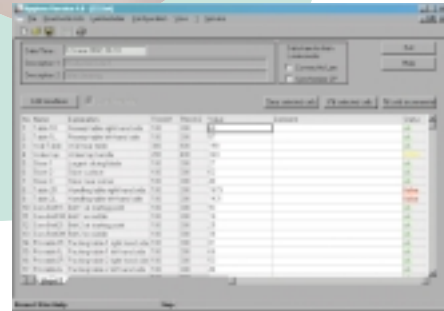

Reagent Options

- **All-in-one reagent devices (Snapshot)**
- **Liquid reagents in dropper bottles**
Contamination Control Kits containing;
Liquid-stable Luciferase, Extractant +/- Pre-wetted extractant swabs

Distributed by



a product of

Hygiena International Limited
Unit 25 Wenta Business Centre
Colne Way Watford
Hertfordshire WD2 4ND UK
Tel: +44 (0)1923 818821
Fax: +44 (0)1923 818825
Email: enquiries@hygiena.net
Web: www.hygiena.net

Hygiena LLC
941 Avenida Acaso
Camarillo
CA 93012-8700 USA
Tel: +1 805 388 8007
Fax: +1 805 388 5531